

# envely

PORTABLE SKIN CARE KIT





For a brilliant young skin  
Total Skin Care Solution – **envely**



# Detailed Names

1 ION/VIBRATION HEAD

2 INFRARED LAMP

3 POWER BUTTON

4 FUNCTIONS

C: CLEANSE

M: MASSAGE

N: NUTRITION

L: LIFTING (Tighten)

W1 : WHITENING (Brighten)

W2 : WRINKLE CARE

5 STEP BUTTON

6 LOW/HIGH BUTTON

7 BATTERY CHARGE LAMP

8 HAND ION CHIP

9 DEVICE HOLDER



## ▪ Main Body Specification

142(H) × 52(W) × 54(T)mm, approx. 102g

Voltage : DC 3.7V, 800mA

## ▪ Adaptor Specification

86(H) × 47(W) × 23(T)mm, approx. 54g

Power : AC 100~240V 50/60Hz

## ▪ Device Holder


157(H) × 67(W) × 34(T)mm, approx. 52g

# How to Use Buttons


You can select function by pressing the buttons.



## ▶FUNCTION BUTTON

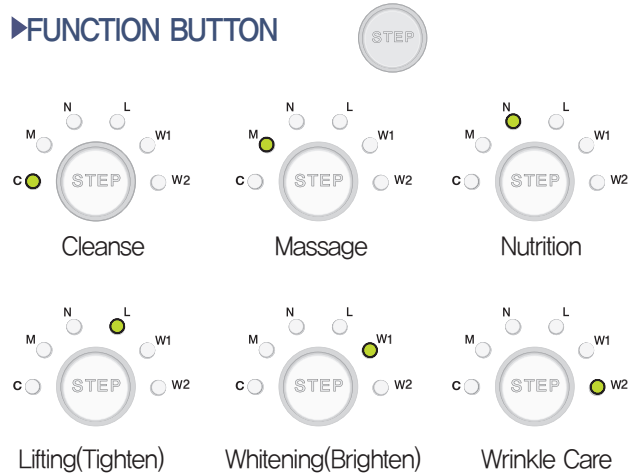

## ▶ L/H button

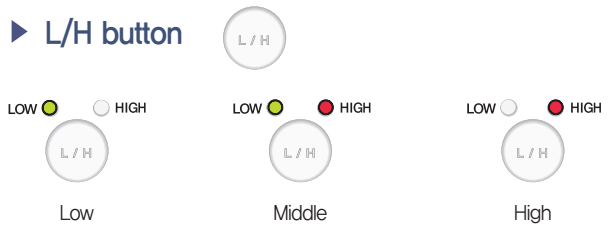


※ Press and hold L/H button to switch off Vibration function.

## Basic Facts about Skin

Your skin is made up of three levels: the **Epidermic layer**, the **Dermic layer** and the **Hypodermic layer**. The Barrier-Zone in the Epidermis protects your skin from the harmful effects of the environment, but can also prevent necessary nutrients from being absorbed.

**envely's** high-end technology uses Galvanic ion, Infrared ray and Vibration at the same time to reduce dead skin cells and sebum. It also increases the absorption of water soluble nutrients and oily nutrients deeper into the skin. Vibration and Infrared ray function helps to improve your skin, making it healthier and more beautiful.

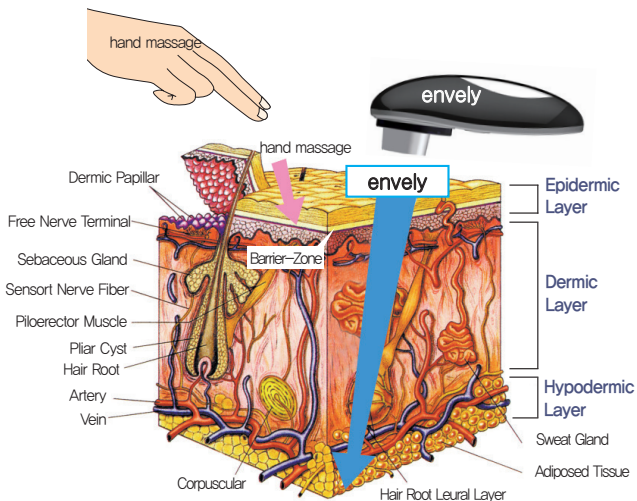

1

### **Galvanic Ion**

Galvanic ion is a safely rectified current. Using both negative and positive ions it cleans and purifies pores and removes dead skin cells. It also transports valuable soluble substances deeper into the skin, keeping your skin young and healthy.

2

### **Infrared Ray**

Infrared ray gives healthy and vibrant skin with its thermal therapy.

3

### **Vibration**

Vibration relaxes and contracts the muscle and helps your skin to be healthier and more elastic.



**3 Functions  
of envely**

# User's Guide

## Normal/Dry Skin


- Rub Continuously  
→ Remain for 3 seconds and move to the next spot



### Step 1. Cleanse (C)



1. Remove makeup and clean your face.
2. Apply cleansing gel, lotion or cream evenly on face.
3. Start evenly by pressing the **POWER** button, then press **STEP** button to select **C(Cleanse)**.
4. Press **L/H** button to select the intensity and hold the device touching ion chip on the back. Next, lay the head on the skin area to be treated and move it softly over the skin in slow circles. (The device is not waterproof)

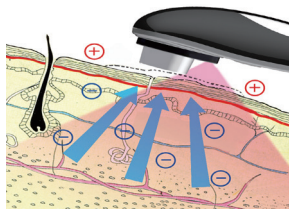


This level ensures cleans and pure skin. The positive ion, infrared ray and vibration function increases positive result. Keratin, sebum and other skin degradation products are simply dissolved and washed away after this treatment.



### Step 2. Massage (M)



1. Apply massage gel, lotion or cream evenly on face.
2. Press **STEP** button to select **M(Massage)**.
3. Hold the device touching ion chip on the back. Next, lay the head on the skin area to be treated and move it softly over the skin in slow circles. Treat once (4minutes) and wipe off with steamed towel or wash face with water.

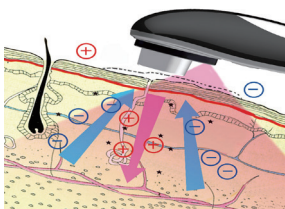


Just like a conventional massage, the evenly massage brings back relaxation and elasticity to the skin. Negative and positive ions, infrared ray and vibration activate the cell spaces and effectively stimulate the metabolism. Your skin beauty improves notably.



### Step 3. Nutrition (N)



1. Apply your daily basis skin care products evenly on the face by hand.
2. Press **STEP** button to select **N(Nutrition)**.
3. Hold the device touching ion chip on the back. Next, lay the head on the skin area to be treated and move it softly over the skin in slow circles.

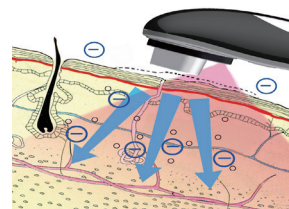

Your skin needs nutrition. The unique functions strongly activates the cell spaces and thus guarantees to transport soluble agents deeper into the skin in a highly efficient way.




TIP 1: After finishing the M(Massage), use skincare products specializing in Anti-wrinkle, Brightening and Tightening for steps N through W2.  
 TIP 2: Use on top of Hydrating or Vitamin C facial mask sheets during N(Nutrition) ~ W2(Wrinkle Care).



#### Step 4. Lifting (L)



1. Do not wipe off the skincare products on your face after **N(Nutrition)**, Press the **STEP button** to select **L(Lifting)**.
2. Hold the device touching ion chip on the back. Next, lay the head on the skin area to be treated and move it softly over the skin in slow circles. Massage in upward direction starting from the chin.

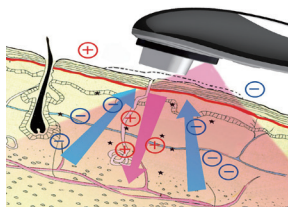

The positive and negative ions with infrared ray and vibration activates the production of elastin and collagen. Tightening and rejuvenation of the skin will make your skin become soft and elastic again.




#### Step 5. Whitening (W1)



1. Apply specialized skin care products (Vitamin C or Brightening ampoules) on face
2. Press **STEP button** to select **W1(Whitening)**.
3. Hold the device touching ion chip on the back. Next, lay the head on the skin area to be treated and move it softly over the skin in slow circles.

※ Concentrated treatment of this step is required on skin with pigments and freckles.

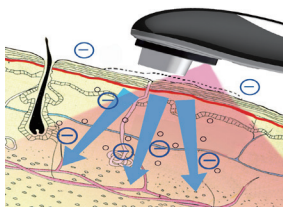


Negative ions, infrared ray and vibration transports essential nutrients for brightening deep into the skin. For a targeted treatment of pigmentary abnormalities and age spots, evenly enables you to experience a notable reduction of skin discolorations.



#### Step 6. Wrinkle Care (W2)



1. Apply anti-age products on the spots to treat.
2. Press **STEP button** to select **W2(Wrinkle Care)**
3. Hold the device touching ion chip on the back. Next, lay the head on the skin area to be treated and move it softly in opposite direction of the wrinkle.

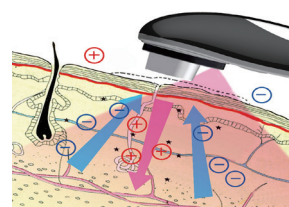

Positive and Negative ions with infrared ray and vibration transports the active components of wrinkle care cosmetic products. By using with special wrinkle care products, the visible effects of tightening and smoothing of wrinkles and sagging skin muscles will be visibly improved.

# User's Guide

## Combination Skin


- Rub Continuously  
→ Remain for 3 seconds and move to the next spot



### Step 1. Cleanse (C)



1. Remove makeup and clean your face.
2. Apply cleansing gel, lotion or cream evenly on face.
3. Start evenly by pressing the **POWER (O)** button, then press **STEP** button to select **C(Cleanse)**.
4. Press **L/H** button to select the intensity and hold the device touching ion chip on the back. Next, lay the head on the skin area to be treated and move it softly over the skin in slow circles. Use this step twice (8minutes).

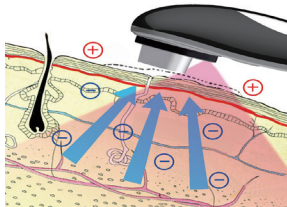


This level ensures cleans and pure skin. The positive ion, infrared ray and vibration function increases positive result. Keratin, sebum and other skin degradation products are simply dissolved and washed away after this treatment.



### Step 2. Massage (M)



1. Apply massage gel, lotion or cream evenly on face.
2. Press **STEP** button to select **M(Massage)**.
3. Hold the device touching ion chip on the back. Next, lay the head on the skin area to be treated and move it softly over the skin in slow circles. Treat once (4minutes) and wipe off with steamed towel or wash face with water.

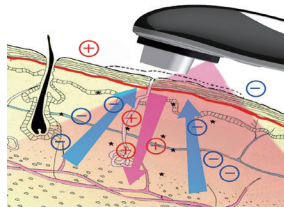


Just like a conventional massage, the evenly massage brings back relaxation and elasticity to the skin. Negative and positive ions, infrared ray and vibration activate the cell spaces and effectively stimulate the metabolism. Your skin beauty improves notably.



### Step 3. Nutrition (N)~ Step 6. Wrinkle Care (W2)



1. Apply hydrating or Vitamin C facial mask sheet on face.
2. Press **STEP** button to select **N(Nutrition)**. Treat over facial mask sheet, (8 minutes)
3. Continue the following steps **L(Lifting)**, **W1(Whitening)** and **W2(Wrinkle Care)**. Peel the the mask sheet off after the care is done.
4. Use your daily basis cosmetic products to finish the treatment.

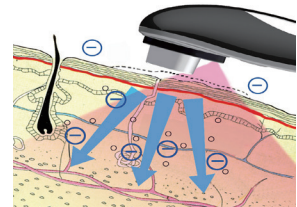

Positive and Negative ions with infrared ray and vibration function transports nutrient essentials in the facial mask sheet deep into the skin and hydrates the skin to become clearer and younger.

# Oily/Troubled Skin


- Rub Continuously
- Remain for 3 seconds and move to the next spot



## Step 1. Cleanse (C)



1. Remove makeup and clean your face as usual.
2. Apply cleansing gel, lotion or cream evenly on face.
3. Start evenly by pressing **POWER** button, then press **STEP** button to select **C(Cleanse)**.
4. Press **L/H** button to select the intensity and hold the device touching ion chip on the back. Next, lay the head on the skin area to be treated and move it softly over the skin in slow circles. Use this step three times. (12minutes) and wipe off with steamed towel or wash face with water.

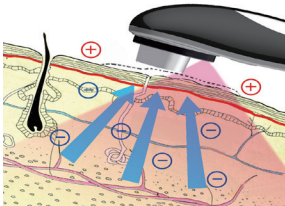


This level ensures clean and pure skin, The positive ion, infrared ray and vibration function increases positive result. Keratin, sebum and other skin degradation products are simply dissolved and washed away after this treatment.



## Step 3. Nutrition (N)~ Step 6. Wrinkle Care (W2)



1. Apply hydrating or Vitamin C facial mask sheet on face.
2. Press **STEP** button to select **N(Nutrition)**. Treat over facial mask sheet twice. (8 minutes)
3. Continue the following steps **L(Lifting)**, **W1(Whitening)** and **W2(Wrinkle Care)**. Peel the mask sheet off after the care is done.
4. Use your daily basis cosmetic products to finish the treatment.

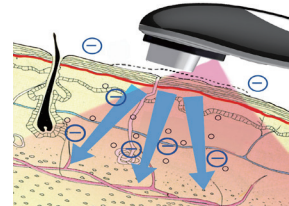

Positive and Negative ions with infrared ray and vibration function transports nutrient essentials in the facial mask sheet deep into the skin and hydrates the skin to become clearer and younger.

# Body Care

## Neck



- Apply neck cream or nutrition, tighten, brighten, wrinkle care product on the neck area.
- Start **envely** by pressing the POWER button(⏻) and press STEP button to select N(nutrition). Press L/H button and select the intensity.
- Hold the device and massage in upward direction, ( ↑↑ )( ↗↗ ). Use Nutrition(N) through Wrinkle care(W2).

## Hand



- Apply massage gel on the back of your hands.
- Start **envely** by pressing POWER button(⏻). Use Cleanse(C) and Massage(M). Use a wet towel to remove the remaining gel.
- Apply hand cream or nutrition cream and select N(Nutrition). Next press the L/H button to select the intensity.
- Use Nutrition(N) through Wrinkle care(W2).

# Please read Before use

---


- | Remove all metallic jewelry before use.
- | **envely** operates only when it contacts with skin.
- | The preset application time for each step is 4minutes. A signal tone sounds after 4 minutes and the device switches to the standby mode. The device will switch-off completely if not reactivated within 2 minutes.
- | If you are treating other person with **envely**, be sure to touch the person with your free hand(make skin contact), in order to close the circuit and ensure proper operation of envely.
- | In case of very sensitive skin or areas around your eyes, gently lay **envely** on the area. Remain for 3 seconds and move to the next spot to treat.
- | For **Dry or Normal Skin**, use envely 3~4 times a week. And for **Sensitive Skin**, use envely 2~3 times a week depending on your skin condition.
- | For **Combination skin or Oily/Troubled skin**, use **envely** 2~3 times a week on the top of Hydrating or Vitamin C facial mask sheets.
- | If you don't experience any irritation you can use it on a daily basis.
- | Use **envely** before bed time for best result of skin regeneration.
- | You might feel slight stinging or irritation when your skin is very dry. The stinging is a normal reaction caused by the ION function. To avoid and reduce irritation, set level to low and apply more gel or cosmetics.
- | Do not use **envely** on children under age of 13.
- | Stop using **envely** if you have the following symptoms and contact your dermatologist;  
Repetitive flushing, itchiness or pyogenic inflammation
- | For allergic or diathesis skin types please talk with your doctor first prior to using this product.

# Maintenance Guide and Recharging envely

- | Do not wash **envely** with water or store in damp places.
- | Clean head with a damp/dry cloth after use.
- | Always keep electrical contacts clean, free of dust and dry.
- | Connect the adaptor to the device after switching–off the power and make sure the **RED** lamp is on. The red lamp turns **GREEN** when the charging process has been completed.
- | If the battery is completely empty, charge for at least 3 hours before using the device again. The device is functional for approximately 5 hours when charged completely.
- | Vibration function operates with the cycle of 3 seconds vibration followed by 1 second pause. Press and hold L/H button to turn off Vibration.
- | The nominal input voltage is 100V–240V AC, meaning that you use the device while travelling everywhere around the world.
- | Remove the adaptor from the device as soon as the device has been charged completely.

## How to check the remained battery

- | Activate **envely** by pressing **POWER** button(⏻), then press and hold the **STEP(function)** button to indicate the charging status of the battery. Device is fully charged when all **STEP(function)** lamps are on.

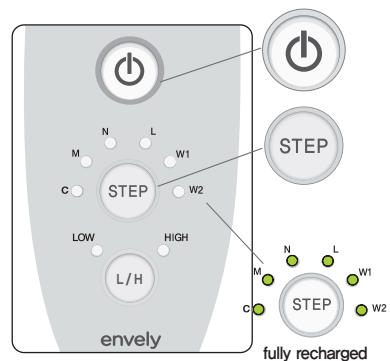



## **envely, a marvel of skin science gets the world's attention**

For a perfect solution and astonishing performance in skin care, Lotts Co.,Ltd. produces **envely** based on years of research and development, **envely** is a skin care device and performs triple functions using galvanic ions, infrared ray and vibration at the same time. **envely** was born with state of the art technology to fulfill the dreams of women all around the world desiring youger and healthier skin.

### **Clinical Tests**

Clinical tests show that when compared to hand massaging alone, the **envely** dramatically increases the number of nutrients and vitamins absorbed into the skin, Vitamins A & C, and Retinol penetrate deep into the epidemic and dermic lasers, discharging sebum, keratin and other harmful waste. The science shows that the **envely** adds nutrition to the skin, firms sagging skin, brightens the face and reduces wrinkles.

### **Skin Care Device Announced to the Academic Journals in US, UK, Japan**

The performance and effects of the clinical tests, were announced in major academic journals. The outstanding effects were reported twice in the SCI Journal, MOLECULAR STRUCTURE 2003 & 2005 in US, six times in Jpan and once in UK. Lotts' tecnology was introduced in a total of nine academic journals specializing in skin care.



envely

## Advanced Technology and Reliability

Lotts' technology is credited and patented around the world.

**envely** is a specialized skin care device with the most advanced technology, reconized by world famous scientific journals for its proven excellence and has acquired various quality marks from domestic and international institutions guarantying quality assurance and customer satisfaction.





# Certifications & Patents



## - Patent -

- Korea Patent No.10-0500100
- Korea Patent No.10-0500911
- Korea Patent No.10-0507453
- Korea Patent No.10-1169675
- Korea Patent No.10-0479622
- China Patent No.ZL 2003 8 0110348,6
- China Patent No.ZL 2003 8 0110349,0
- Japan Patent No.43519593
- Japan Patent No.4431536
- US Patent No.US 7,922,675 B2
- US Patent No. US 8,591,437 B2
- US Patent No. US 8,585,619 B
- Europe Patent No. 2471453

## - Clinical Test -

- Ion Induction, Ion Deduction, Dual Effect of Ion & Ultrasonic Wave performed by Hiroshima Prefectural University, Japan
- Ion Induction, Ion Deduction, Dual Effect of Ion & Ultrasonic Wave performed by Mayssei University, Japan
- In October, 2008, Kangdong Sacred Heart Hospital
- In November 2008, Guri Hospital of Hanyang University Semyong University, Korea

## - Jounal/Thesis/Books -

- 2002 FRAGRANCE JOURNAL (Japan, 2 times)
- 2003 FRAGRANCE JOURNAL (Japan)
- 2003 MOLECULAR STRUCTURE JOURNAL (US)
- 2003 ICAVS-2 (UK)
- 2003 Skin Protection and Biotechnology (Japan)
- 2004 How to Become Young-Reality (Japan)
- 2005 MOLECULAR STRUCTURE JOURNAL (US)
- 2005 Beauty and Skin Science Journal (Japan)

## - Certification-

- KC (MSIP-REM-LMe-ENVELY)  
CE (EN55014, EN61000)
- FCC (ANSI C63.4 :2009)
- RoHS (IEC 62321, EN 50581)
- PL Insurance (Korea Federation of Small and Medium Business)
- ISO 9001 : 2008
- Korean Testing Certification (ion, vibration, infrared ray)



KC(Korea Certification)Mark



Mark for Hygiene and  
Quality Hygiene



CEMark



FCCMark

**RoHS**

RoHS Mark



ISO Mark

The leader in portable skin care kit



**Lotts**

---

Lotts Co.,Ltd. #1005, 10F, LG Twintel, 6, Samseong-ro 96 gil, Gangnam-gu, Seoul, Korea  
Tel : 02-565-7500(代) Fax: 02-561-1472 Customer: 02-565-8787  
Factory : 69, Techno 8-ro, Yuseong-gu, Daejeon, Korea  
www.lotts.co.kr E-mail: lotts@lotts.co.kr